

EUDRACAP® colon

From the upper small intestine to the ileo-colonic region.
Optimize disease treatment efficiency with targeted drug delivery.

Targeted drug delivery to the ileo-colonic region

EUDRACAP® colon is designed for drug products targeted for release in the ileo-colonic region. Like EUDRACAP® enteric, it provides reliable gastric resistance to protect sensitive ingredients.

Reliable protection and delayed release


Dissolution testing of EUDRACAP® colon confirms that no drug release occurs during the acid stage (HCl 0.1N) or at buffer stage at pH 6.8. After conversion to the buffer stage at pH 7.2, release starts.

In comparison to EUDRACAP® enteric, which shows fast release at pH 6.8, EUDRACAP® colon provides significantly more delayed release at a higher pH.



For sensitive drug actives

EUDRACAP® colon provides reliable protection and targeted delivery of sensitive actives like live biotherapeutics, polypeptides and proteins (biologics), and nucleotides.

Target release with EUDRACAP® colon and EUDRACAP® enteric

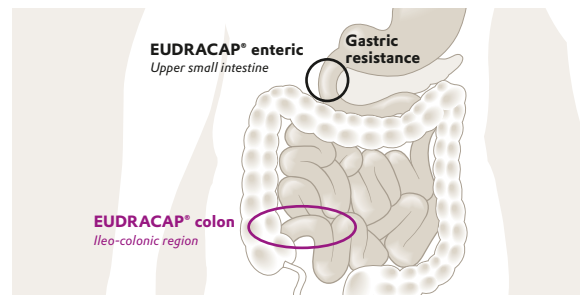

Dissolution profile using caffeine as model compound

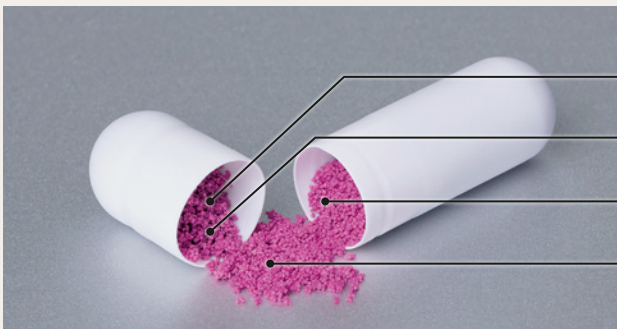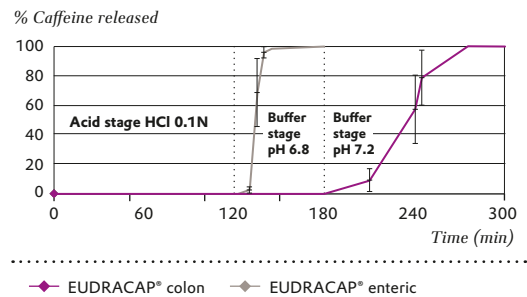

- Live biotherapeutics
- Biologics/proteins
- Nucleotides
- Other sensitive actives

This information and all further technical advice are based on our present knowledge and experience. However, it implies no liability or other legal responsibility on our part, including with regard to existing third party intellectual property rights, especially patent rights. In particular, no warranty, whether express or implied, or guarantee of product properties in the legal sense is intended or implied. We reserve the right to make any changes according to technological progress or further developments. The customer is not released from the obligation to conduct careful inspection and testing of incoming goods. Performance of the product described herein should be verified by testing, which should be carried out only by qualified experts in the sole responsibility of a customer. Reference to trade names used by other companies is neither a recommendation, nor does it imply that similar products could not be used.

Evonik Operations GmbH
Health Care Business Line

healthcare@evonik.com
evonik.com/healthcare